



Comfort Series 500

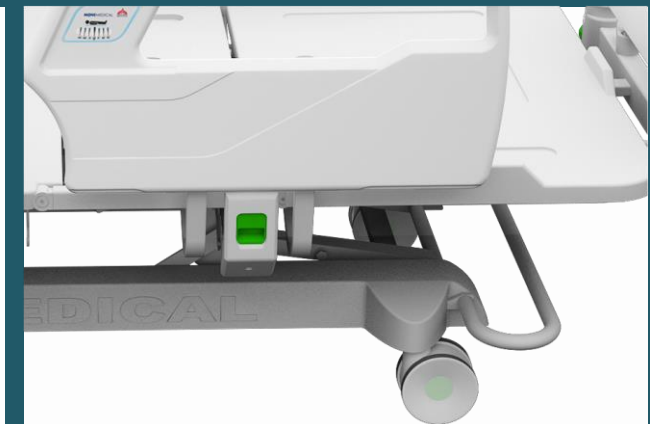

Medispecs Comfort Series 500 Incorporates 4 electronic motors that work together to control height, tilt, back rest, knee-break and calf section adjustment.

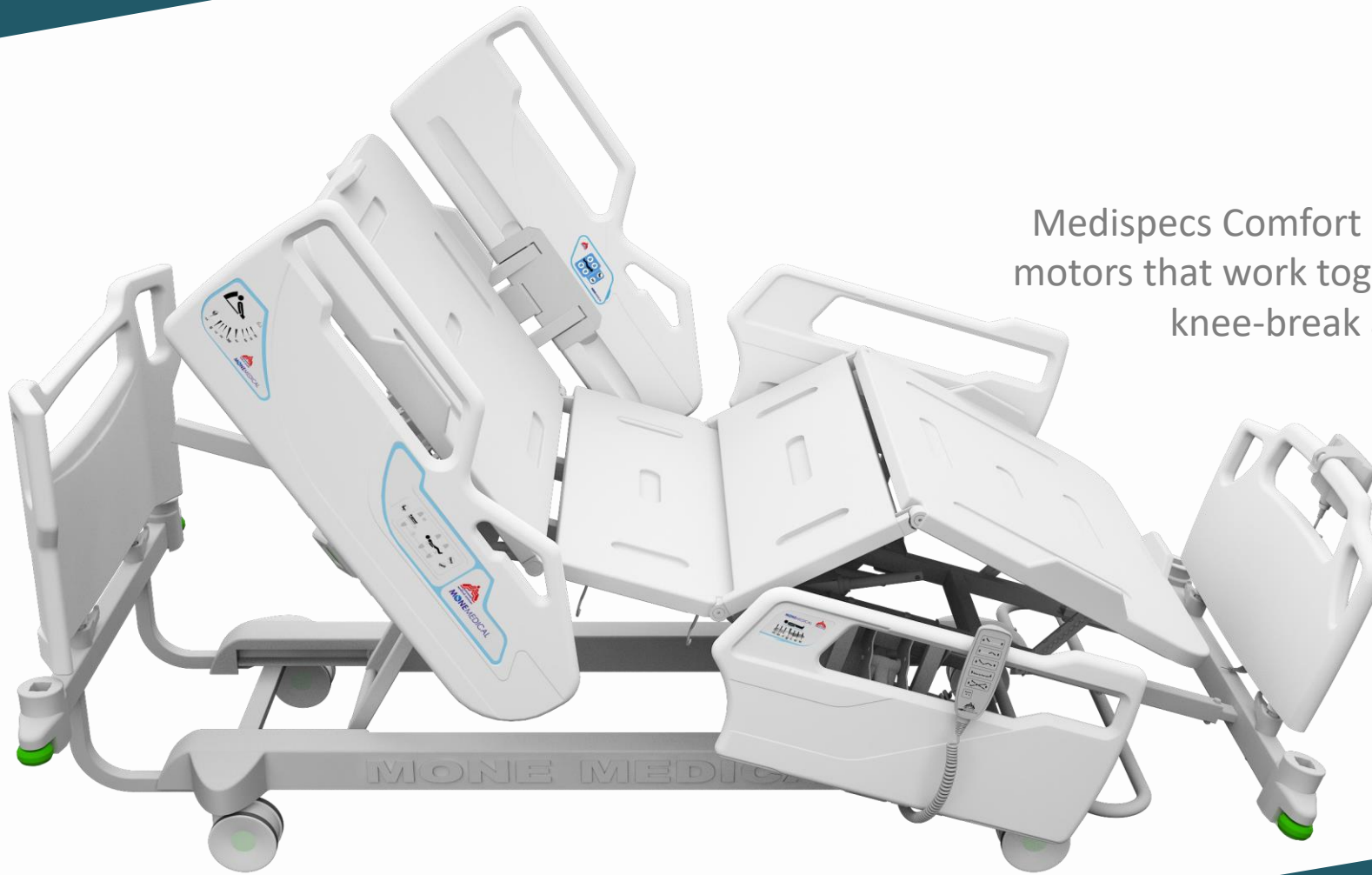

Comfort Series 500 HB06-S

Medispecs Comfort series 500 has been designed for use in ICU and High Care environments. This bed has been designed to deliver both caregiver and patient optimal ease of use and comfort with the priority on safety. Offering both handheld adjustment as well as caregiver side panel controls along with underbed lighting helps exceed stringent quality and safety standards while providing superior care.

Key Features:

- Electrically operated height, tilt, back rest, knee-break and calf section adjustment with electrically operated Trendelenburg and reverse Trendelenburg.
- Electronic Auto CPR and manual CPR Function.
- Built in weight scale accurate up to 0.5kg.
- Built in bed exit alarm.
- Underbed night light.
- Built in X-Ray tray in bed head section.
- Antimicrobial, waterproof mattress.
- 4 hook Intravenous pole.
- Angle Indicator for back rest and thigh section.
- Side rails are removable for easy cleaning and manufactured from vacuum formed high impact plastic.
- Bed control functions are incorporated within the side rails for both patient and Carer.
- Illuminated hand remote control for both patient and Carer.
- Full function lockouts for Carer.
- Calf section can be lowered and raised manually by Carer.
- Safe working load of 250kg.
- Variable height of 37cm to 77cm.
- Anti Static directional castors with linked braking system.
- Customisable trim plates. (Colours/Logos/Patterns)

Technical Specifications

Width	: 104 cm
Length	: 220 cm
Width (Lying Surface)	: 90 cm
Length (Lying Surface)	: 196 cm
Back Segment Length	: 68 cm
Back Segment Angle	: 0-70°
Height Adjustment	: 37-77 cm
Foot Segment Angle	: 0-45°
Trendelenburg	: 17°
Reverse Trendelenburg	: 17
Safe Maximum Load	: 250 kg
Shipment Dimension	: 220 x 113 x 90 cm



Standard Accessories

- Infusion Pole
- Drainage Pole
- Mattress
- Manual CPR
- Bed Exit alarm
- Weight Scale
- X ray Compartment

Guarantee



Quality

