

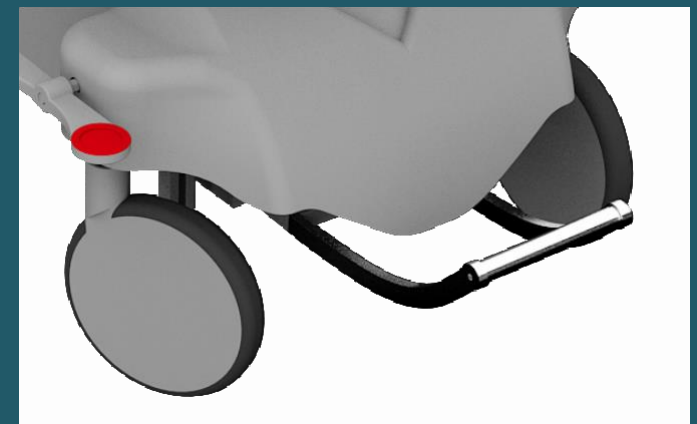


# Emergency Response Table

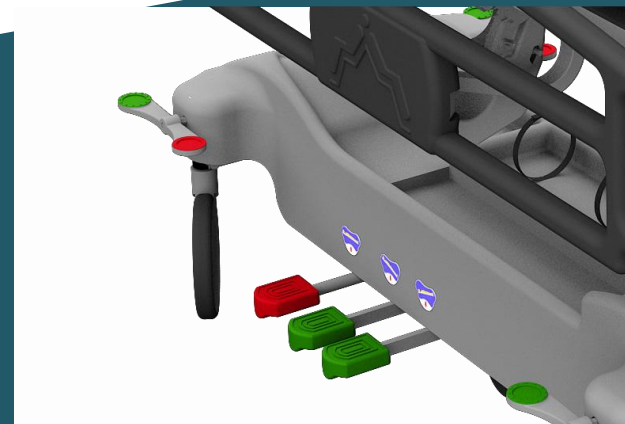

Medispecs Emergency Response Table incorporates manual hydraulic actuators that work together to control height, tilt, back rest and Trendelenburg adjustments

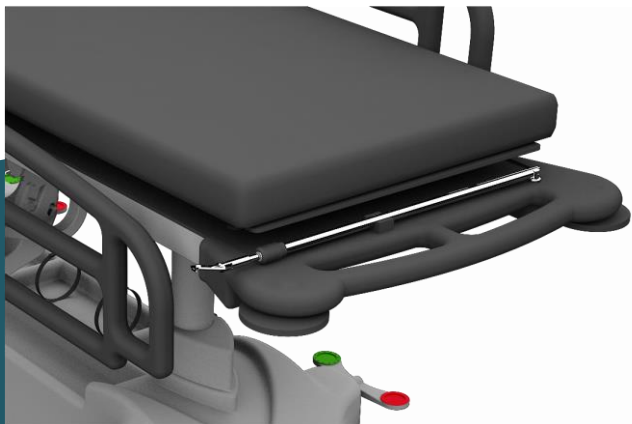

## Emergency Response Bed HS01-S

Medispecs multipurpose hydraulic emergency response table is used for emergency and recovery, patient transportation and minor procedures. The safe and sturdy construction makes this bed ideal for heavy use. The versatile range of options and adjustments ensure efficient patient transfer and ergonomic working conditions for health care staff.


### Key Features:

- Hydraulically operated height, tilt, back rest, Trendelenburg and reverse Trendelenburg.
- Manual adjustable back rest up to 80 degrees.
- Built in X-Ray imaging platform.
- Antimicrobial, seamless waterproof mattress.
- 4 hook Intravenous pole..
- Collapsible side rails are removable for easy cleaning and manufactured from vacuum formed high impact plastic.
- Earthed socket connections for defibrillator and other electronic devices.
- Secure oxygen storage on the sub frame.
- Safe working load of 250kg.
- Variable height of 65cm to 93cm.
- Large Anti Static directional castors with independent and linked braking system.

### Technical Specifications

<b>Outside Dimensions</b>	
Width	: 75 cm
Length	: 220 cm
<b>Lying Surface</b>	
Width	: 62 cm
Length	: 190 cm
Back Segment Angle	: 0-80°
Height Adjustment	: 65-93 cm
Trendelenburg	: 17°
Reverse Trendelenburg	: 17°
Safe Maximum Load	: 250 kg

### Positions of The Bed

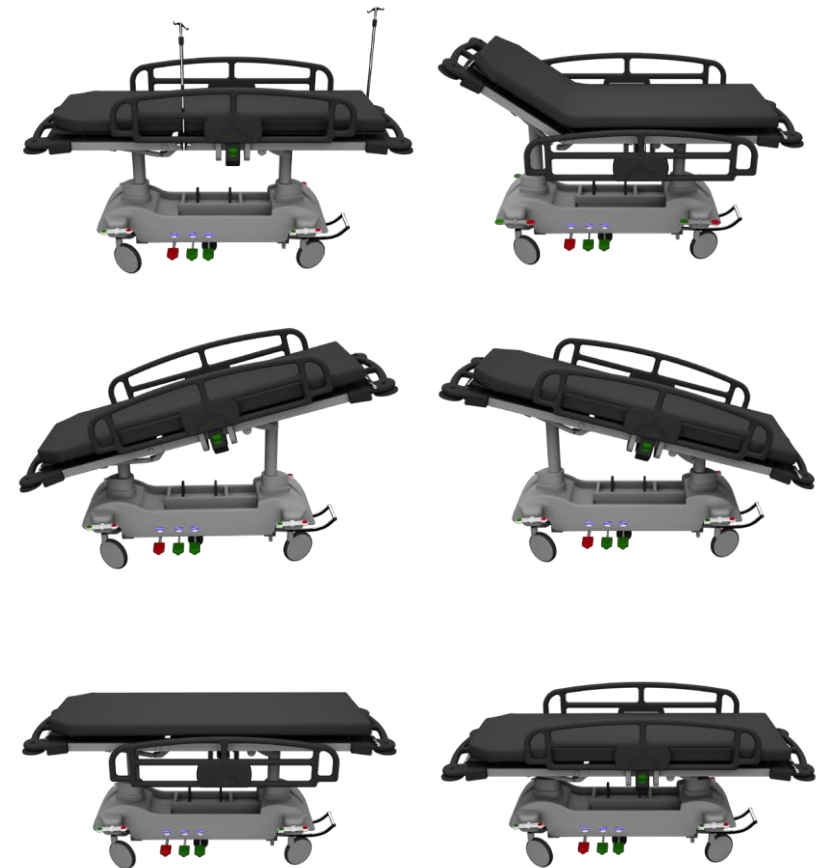


### Standard Accessories

- IV Pole
- Mattress

### Optional Accessory

- 5<sup>th</sup> Castor for manoeuvrability.

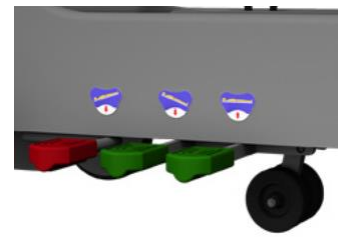

### Technical Support



### Quality

