




## MONESTAR DIALYSIS & CHEMOTHERAPY CHAIR DA07-FW

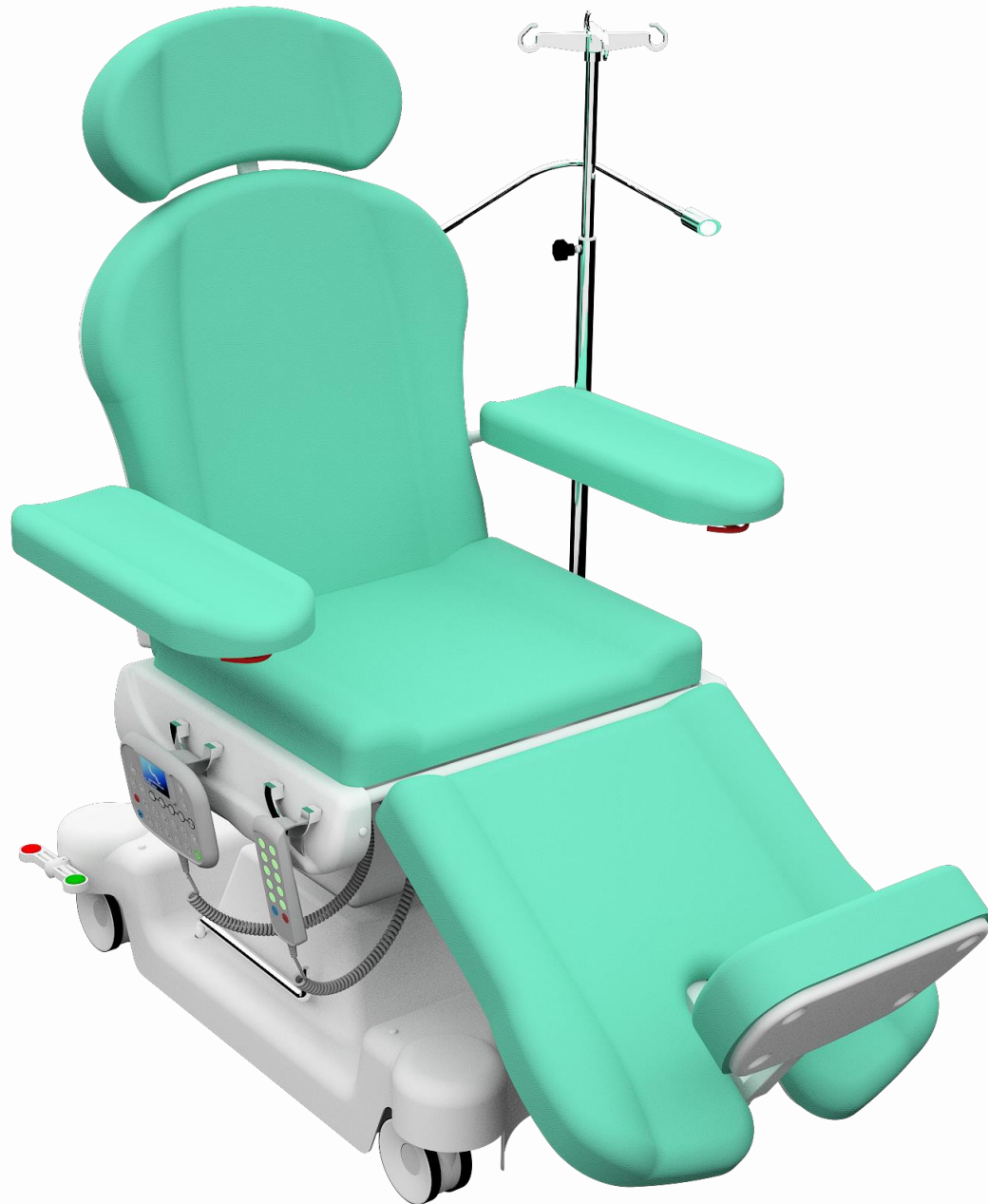

Medispecs Monestar Dialysis chair combines functionality with comfort for both the patient and caregiver. This chair is user friendly and exceptional in its engineering capabilities and functionality. The illuminated hand control, arm rests and five electric motors allow for quick positioning to serve all needs and purposes.

### Key Features

- Electrically operated height, tilt, back rest, knee-break and foot section adjustment with electrically operated Trendelenburg and reverse Trendelenburg.
- Electronic Auto CPR/Emergency function.
- Built in electronic weight scale.
- Antimicrobial viscoelastic foam is used for the padding. This seamless padding moulds itself to the body providing additional comfort and support for long periods. Seamless padding reduces the risk of cross contamination by eliminating sewing and creating a smooth clear surface area which is easy to clean and disinfect.
- 4 hook Intravenous pole
- Adjustable reading lamp
- Illuminated hand remote control for both patient and Carer.
- Full function lockouts for Caregiver.
- The detachable foot rest and head rest is fully adjustable to accommodate the height of the patient.
- Safe working load of 230kg.
- Variable height of 55cm to 75cm.
- Anti Static directional castors with independent braking system.
- This multi purpose chair can be used in numerous disciplines and procedures such as Dialysis, Chemotherapy, radial catheterization, Oncology and Phlebotomy



**Technical Specifications**

Width	: 110 cm
Length	: 210 cm
Back Segment Length	: 67,6 cm
Back Segment Angle	: 0-75°
Min-Max	: 55-75 cm
Seat Length	: 51,7 cm
Foot Segment Length	: 62,5 cm
Foot Cramp Support	: 20 cm
Foot Segment Angle	: 0-42°
Trendelenburg	: 0-14°
Safe Maximum Load	: 230 kg
Net Weight Gross	: 114 kg
Weight Shipment	: 174 kg
Dimension	: 220x113x90 cm

**Standard Accessories**

- Nurse Control Panel
- Electronic weight scale
- Infusion Pole
- Headrest
- Reading Lamp

**Optional Accessories**

- Cushions
- Towel-Holder
- Back up Battery
- LCD Screen
- Dining Table

**Color Options**



**Optional Accessory**



**Technical Support**



**Quality**

